

MEDSTOCK[®]

**TRANSPARENT
FILM**

DRESSING

**A HIGHLY PERMEABLE
POLYURETHANE FILM DRESSING**

Medstock Transparent Film Dressing is a highly permeable polyurethane (PU) film dressing. It is made of a medical-grade, pressure-sensitive adhesive coating.

AVAILABLE IN 4 SIZES

Key Functions and Features


- + Waterproof, flexible and conforms well to body contours to enhance patient comfort
- + Prevents bacterial invasion and accelerates wound healing
- + Provides protection against friction and shearing
- + Highly breathable

General Advice on Safe Use

- + Ensure wound area is clean and dry before application
- + Can be left in place for up to 7 days depending on exudate level and type of wound
- + Can be used as a secondary dressing to secure other primary dressings

Indications for Use

- + Suitable for wounds with nil to low exudate
Ideal for epithelialising wounds
- + Waterproof film can be used to secure the fixation of a wide variety of dressings, tubes, catheters, drains and other medical devices
- + Can be used directly over superficial wounds (Superficial burns/scolds, abrasions, grazes and cuts)
- + Stage I and II pressure injuries
- + Can be used over surgical wounds to spread force and reduce suture load

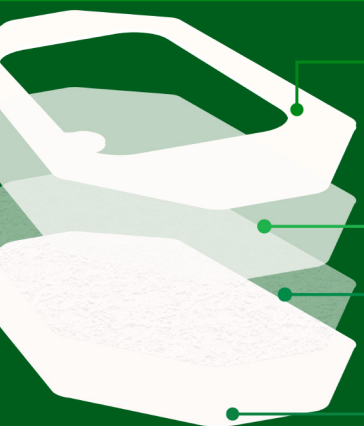

First Layer: CCK paper

Second Layer: Polyurethane (PU) film

Third Layer: Pressure sensitive adhesive

Backing Release Paper: CCK paper



Contraindications / Precautions

- + Seek medical advice immediately for all serious wounds and burns
- + If reddening or sensitisation occurs, discontinue use and consult a physician

Storage Conditions

- + Store in a well-ventilated clean room at or below 35°C with relative humidity less than 80% and no corrosive gas

Period of Validity

The product shall be valid for 5 years from the sterilisation date, provided that it is kept under the stipulated storage conditions

Packaging

- + Each dressing is packaged in a paper/poly closed pouch, then packed into our unique colour coded Medstock boxes

Care and Use Symbols



Read Instructions for use.



Single use only.



Product is latex free.



Sterilised using ethylene oxide.



Keep dry.



Store out of direct sunlight.



Do not use if packaging is broken or damaged.

Product Code	Size	Pcs/Box
MS100607TF	6cm x 7cm	100pcs/box
MS1001012TF	10cm x 12cm	50pcs/box
MS1001520TF	15cm x 20cm	50pcs/box
MS10010130TF	10cm x 30cm	50pcs/box